

More Control in Diabetes Management.




GlucoBalance™ Data Management Software

GlucoBalance can help you or your patient manage diabetes more closely by tracking and analyzing blood glucose results. It enables transfer of blood glucose results stored in the memory of the meter directly to a personal computer, enabling review and analysis of test results using a variety of charts and graphs.

- Encourages tighter control in diabetes management
- Electronic log book eliminates the need to record test results by hand
- Useful in tracking the number of tests performed daily and recognizing trends
- Visual tools provide an easy-to-understand summary of blood glucose results

GlucoBalance may be used in conjunction with the following Hypoguard blood glucose meters:

- QuickTek™
- Hypoguard Advance®

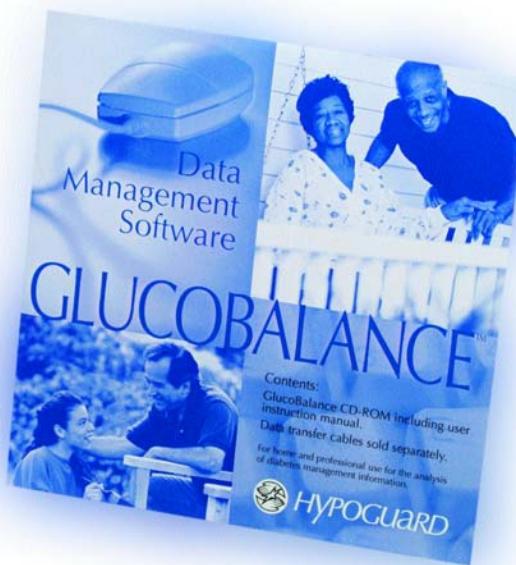

GLUCOBALANCE™

A Valuable Tool.

You can provide a valuable service to your patients by meeting with them regularly to review their results. A variety of chart formats are available to meet your needs.





Electronic Log Book Report

Makes it easy for healthcare professionals to view a patient's history in a diary-like format.

Pie Chart Report

Each "slice" of the pie is color-coded to represent a blood glucose range.

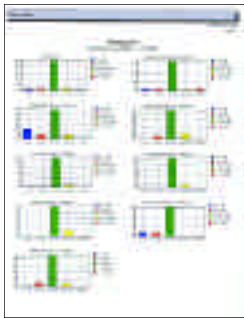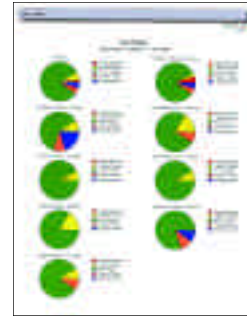

Histogram Report

A color-coded bar chart format for different times of the day. Presents data in separate views of each timeslot and in an aggregate view of all timeslots.




Statistics Report

Provides the number of high readings, low readings, the mean, and standard deviation for the time period selected.


Modal Day Report

Readings and averages are shown across each timeslot for the chosen time period. Target blood glucose levels are represented in green. Hyperglycemic and hypoglycemic levels are displayed with red and blue lines.

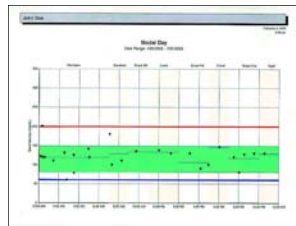

Test Frequency Report

Indicates all readings by timeslot for each day in the defined time period. It computes the average number of tests per day, the actual number of readings per day, and the percent of tests taken within each timeslot.

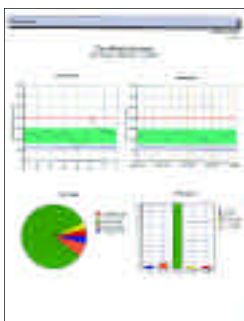


Two Week Summary Report

Combines two weeks worth of information showing pie and trend charts, the modal day and histogram reports.

Trend Chart Report

Makes it easy to compare actual tests against the target zone to trend daily fluctuations or triggers.

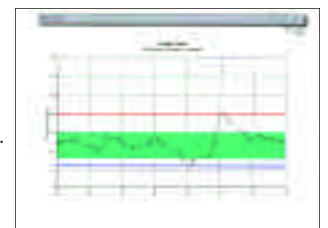

Call 1-800-566-8558 to order GlucoBalance™ software and meter cables.

Product Name	Product Number
GlucoBalance Data Management Software	555030
GlucoBalance Download Cables for Hypoguard Meters	
Cable for QuickTek™	555020
Cable for Hypoguard Advance®	555022

Patients should always remember to consult with health care professionals before making changes to their diabetes management routines.



Minneapolis, MN 55439 USA
 Customer Service: 800-566-8558
 Fax: 952-646-3110
www.hypoguard.com